BENEFIT[®]

Leading in Orthodontic Implants

psm-medical.com

Quality is never an accident.

It is always the result of high intention, sincere effort, intelligent direction and skillful execution; it represents the wise choice of many alternatives.

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Introduction | Einführung

Skeletal anchorage in orthodontics

An important requirement for successful orthodontic treatment is stable anchorage.

Conventional orthodontic anchorage can be limited by periodontal disease and/or missing teeth. Skeletal anchorage has been proven to be very effective. It can be an essential part of non compliance therapy and broadens the options available to the orthodontic practitioner.

Mini Implants are especially advantageous because of the minimally invasive insertion protocol and their relatively low cost.

However, limitations can occur when conventional Mini Implant Systems are used for orthodontic anchorage. In collaboration with leading specialists, **psm** has developed the unique **BENEFIT®** TAD System for temporary fixation in the palate. **BENEFIT®** offers high precision titanium implant screws and a large variety of pre-fabricated connecting parts such as **BENEplates®**, abutments and palatal expanders.

With the introduction of the **BENEFIT®** System **psm** has considerably enhanced the therapeutic diversity of the treatments available for Mesialization, Distalization, RPE as well as uprighting or anchorage of the dentition.

With **BENEFIT®**, **psm** was the first company to offer a full range of skeletal anchorage appliances for the palate – a trend that has become „state of the art“ in Orthodontics.

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Same instruments can be used for several **psm** products. In order to find the right instruments for your use, they are indicated in colour.

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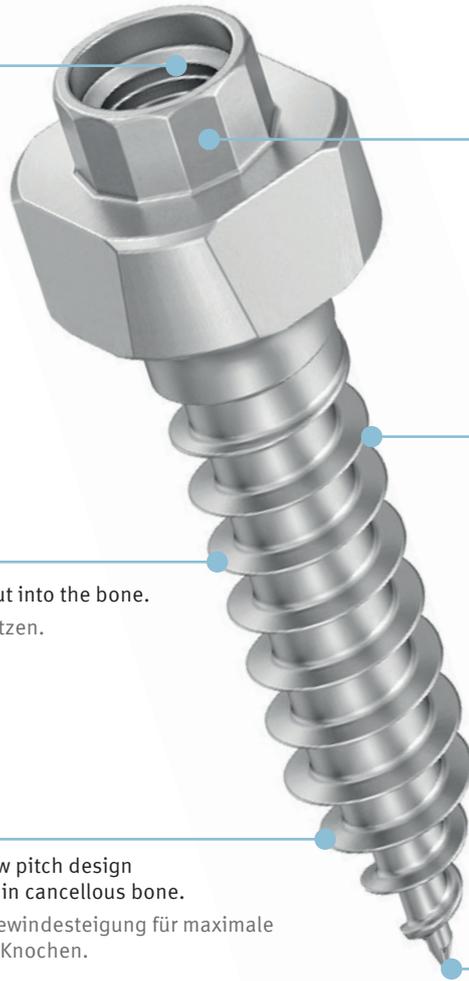
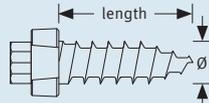
BENEFIT® Orthodontic Screw 1.5mm, 1ea.
BENEFIT® KFO-Schraube 1,5mm, 1Stk.

 All **BENEFIT®** screws now also available as a pack of 2

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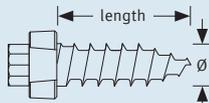
BENEFIT® Orthodontic Screw 1.8mm, 1ea.
BENEFIT® KFO-Schraube 1,8mm, 1Stk.


2x
ST-33-54XXX-2
BENEFIT® Screws in a package of 2ea.
BENEFIT® Schrauben in VPE 2Stk.

Ø2.0

 Size (Ø x length)
 Abmessung (Ø x Länge)

ST

ST-33-54205	2.0 x 5mm	NEW
ST-33-54207	2.0 x 7mm	
ST-33-54208	2.0 x 8mm	NEW
ST-33-54209	2.0 x 9mm	
ST-33-54210	2.0 x 10mm	NEW
ST-33-54211	2.0 x 11mm	
ST-33-54213	2.0 x 13mm	
ST-33-54215	2.0 x 15mm	

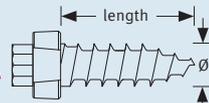BENEFIT® Orthodontic Screw 2.0mm, 1ea.
BENEFIT® KFO-Schraube 2,0mm, 1Stk.

10x

ST-33-54207-10	2.0 x 7mm
ST-33-54208-10	2.0 x 8mm
ST-33-54209-10	2.0 x 9mm
ST-33-54210-10	2.0 x 10mm
ST-33-54211-10	2.0 x 11mm


BENEFIT® Screw dispenser with 10ea.
BENEFIT® Schrauben im Spender 10Stk.

Ø2.3

 Size (Ø x length)
 Abmessung (Ø x Länge)

ST

ST-33-54307	2.3 x 7mm
ST-33-54309	2.3 x 9mm
ST-33-54311	2.3 x 11mm
ST-33-54313	2.3 x 13mm
ST-33-54315	2.3 x 15mm

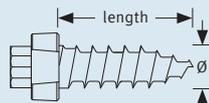BENEFIT® Orthodontic Screw 2.0mm, 1ea.
BENEFIT® KFO-Schraube 2,0mm, 1Stk.

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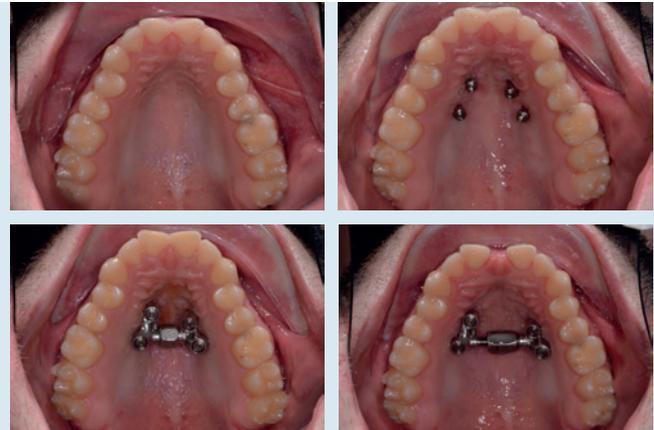
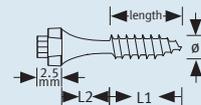
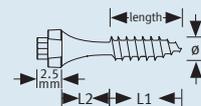
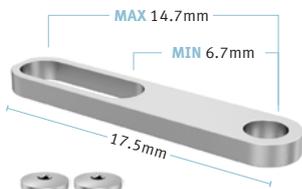
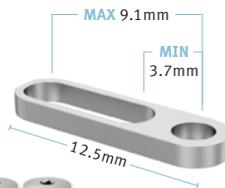
BENEplate | BENEplate

33-54472 1 ea. / 1Stk.

BENEplate, L, paramedian, incl. 2 fixation screws, with 1.1mm wire (12cm), stainless steel

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33-54471 1 ea. / 1Stk.

BENEplate, S, paramedian, incl. 2 fixation screws, with 1.1mm wire (12cm), stainless steel

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33-54422 1 ea. / 1Stk.

BENEplate, XS, paramedian, incl. 2 fixation screws, with 1.1mm wire (12cm), stainless steel

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1.1mm wire
1,1mm Draht

33-54407 1 ea. / 1Stk.

BENEplate, L, incl. 2 fixation screws, with bracket, stainless steel

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33-54474 1 ea. / 1Stk.

BENEplate, S, paramedian, incl. 2 fixation screws, with 0.8mm wire (12cm), stainless steel

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33-54480


1ea. / 1Stk.

BENEplate, L, median, incl. 2 fixation screws, with central 1.1mm wire (12cm), stainless steel, acc. Dr. Nienkemper, Düsseldorf, Germany

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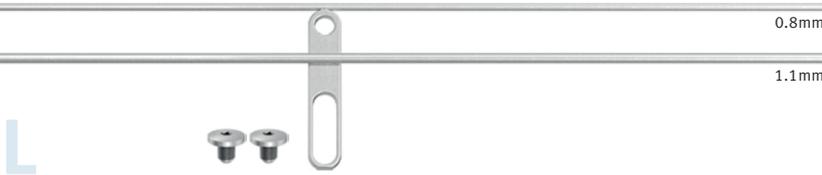33-54485


1ea. / 1Stk.

BENEplate, L, T-Mesialslider plate, median, incl. 2 fixation screws, with central 1.1mm wire + 0.8mm wire anteriorly (12cm), stainless steel

BENEplate, L, T-Mesialslider plate, median, inkl. 2 Fixierschrauben, mit 1,1mm mittigem Draht + 0,8mm Draht anterior (12cm), Implantatstahl

33-54485

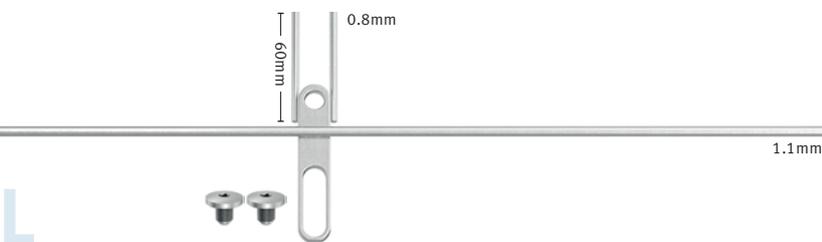33-54490


1ea. / 1Stk.

BENEplate, L, T-Mesialslider plate, median, incl. 2 fixation screws, with central 1.1mm wire (12cm) + 0.8mm wire laterally (6cm/ea.), stainless steel

BENEplate, L, T-Mesialslider plate, median, inkl. 2 Fixierschrauben, mit 1,1mm mittigem Draht (12cm) + 0,8mm Draht lateral (je 6cm), Implantatstahl

33-54493

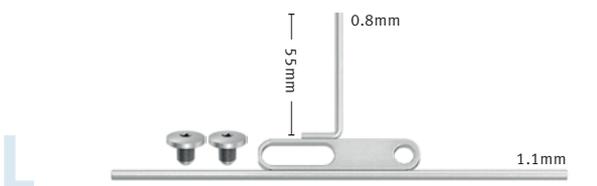
33-54492

1ea. / 1Stk.

BENEplate, L, T-Mesialslider plate, paramedian, incl. 2 fixation screws, with 1.1mm lateral wire (12cm) + 0.8mm wire vertically, stainless steel

BENEplate, L, T-Mesialslider plate, paramedian, inkl. 2 Fixierschrauben, mit 1,1mm Draht (12cm) lateral + 0,8mm Draht vertikal, Implantatstahl

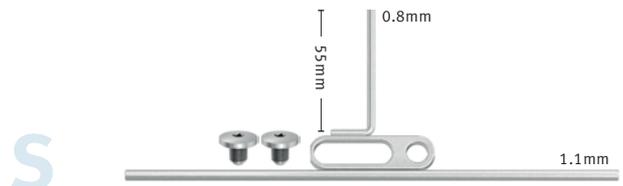33-54493

1ea. / 1Stk.

BENEplate, S, T-Mesialslider plate, paramedian, incl. 2 fixation screws, with 1.1mm lateral wire (12cm) + 0.8mm wire vertically, stainless steel

BENEplate, S, T-Mesialslider plate, paramedian, inkl. 2 Fixierschrauben, mit 1,1mm Draht (12cm) lateral + 0,8mm Draht vertikal, Implantatstahl

BENEplate | BENEplate

Video


33-54730

1ea. / 1Stk.

Arnoudizer for indirect anchorage, incl. 1 fixation screw, stainless steel, acc. Dr. A. Schreurs, Genk, Belgium

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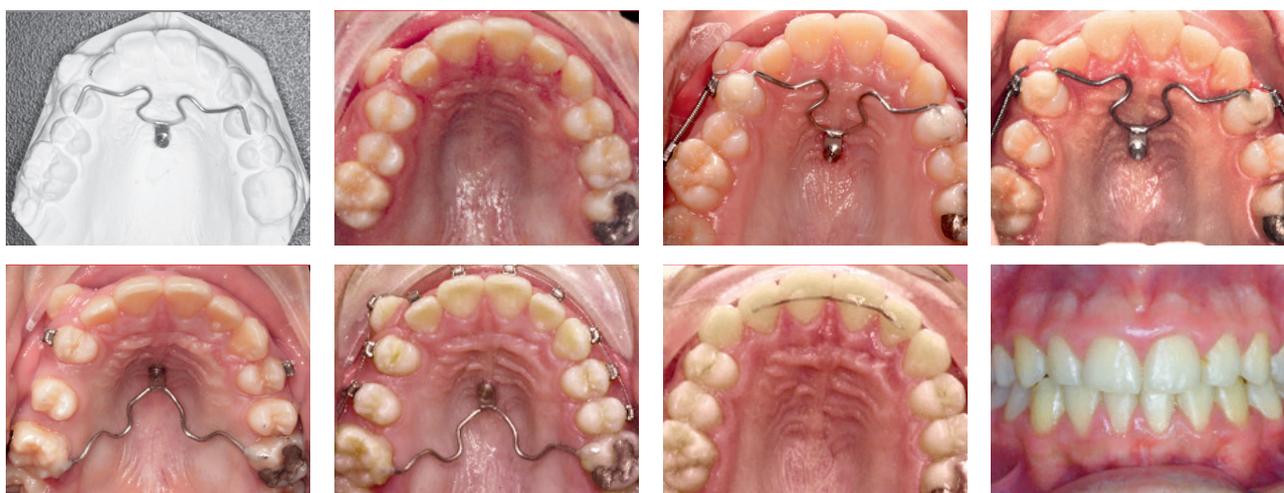A device for molar distalization. The upper molars will be derotated, uprighted and distalized in only two or three months. By converting the system in a transpalatal bar, it can be used as an anchor for the molars. The anterior teeth can be retracted and aligned simply.

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33-54465

1ea. / 1Stk.

Hyrax® Ring, incl. 1 fixation screw, with 1.1mm wire (6cm) laterally, stainless steel

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
33-54461

1ea. / 1Stk.

Hyrax® Ring, incl. 1 fixation screw, with 1.1mm wire (12cm), stainless steel

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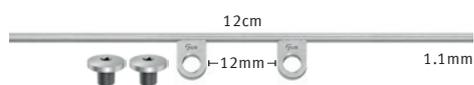
33-54464

1ea. / 1Stk.

Hyrax® Ring, incl. 2 fixation screws, with 1.1mm wire (12cm), stainless steel

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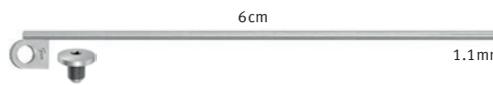33-54750

1ea. / 1Stk.

Hyrax® Ring, incl. 1 fixation screw, with 1.1mm wire (6cm), stainless steel

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33-54535

2ea. / 2Stk.

BENEsTube, standard
BENEsTube, standard

33-54524 240 gr.

33-54525 500 gr.

BENEslieder, NiTi springs, 2 ea.
BENEslieder, NiTi Federn, 2 Stk.

33-54536

2ea. / 2Stk.

BENEsTube, standard
BENEsTube, standard

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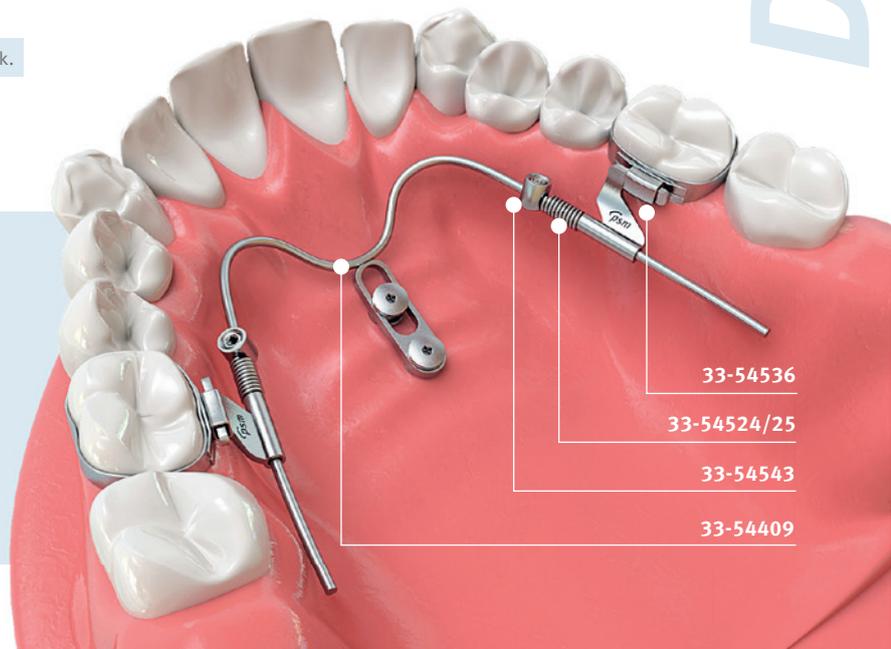
The **BENEs**lider

The most classical indication for skeletal anchorage using Mini-Implants with abutments is molar-distalization and/ mesialization in the maxilla. By means of the so called “**BENEs**lider” the molars can be bodily distalized and/or mesialized using the Mesialslider. In many cases, tooth extraction or dental implants can be avoided.

To couple two **BENEFIT**® Mini-Implants, a **BENEplate** with prewelded wire is connected to the implants and the **BENEs**lider parts (mobilizer, springs, **BENEtube**) are placed on the wire. Active force is applied by springs (240g or 500g) and activated with the Mobilizer.

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BENEFIT® Abutment with 1 bracket
BENEFIT® Abutment mit 1 Bracket

33-54430 1ea. / 1Stk.

BENEFIT® Abutment Standard
BENEFIT® Abutment Standard

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

33-54466 1ea. / 1Stk.

BENEFIT® PEEK Abutment, incl. 1 fixation screw
BENEFIT® PEEK Abutment, inkl. 1 Fixierschraube

33-54469 1ea. / 1Stk.

BENEFIT® PEEK Abutment, 15° angled, incl. 1 fixation screw
BENEFIT® PEEK Abutment, 15° abgewinkelt, inkl. 1 Fixierschraube

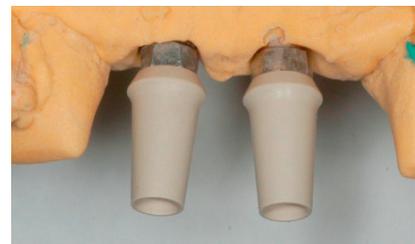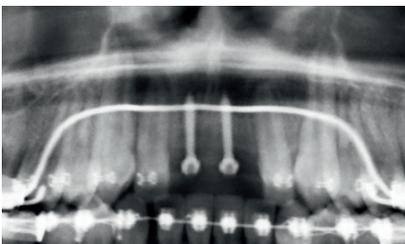Hybrid Hyrax® RPE* | Hybrid Hyrax® GNE*


The **Hyrax®** Ring and **Hyrax®** Abutment facilitate the construction and fixation of a **Hybrid Hyrax®** appliance. It allows welding or soldering with the expansion screw as well as little rotational movements when connecting the appliance to the **BENEFIT®** Screws.

The **Hyrax®** Ring is supplied with 1 fixation screw each.

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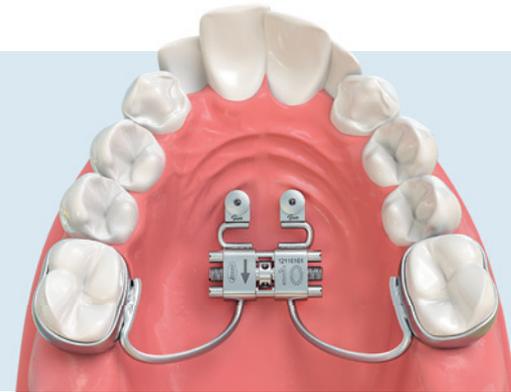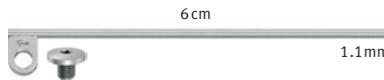
33-54463 1ea. / 1Stk.

BENEFIT® Hyrax®, Abutment
BENEFIT® Hyrax®, Abutment

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bmX® – BENEFIT® Maxillary Xpander | bmX® – BENEFIT® Maxillary Xpander

After using the **BENEFIT®** Screw for the **Hybrid Hyrax®** appliance very successfully, in the last years the **bmX®** now follows the Orthodontists request for an entirely **BENEFIT®** mini-implant anchored appliance for rapid maxillary expansion on young patients. Using the De Franco parallel guide, (Page 35) the distance and the parallelity of the **BENEFIT®** Screws can be exactly determined, on the patient directly. An impression is no longer necessary. Alternatively, visualisation can, of course, be made by means of digital imaging or traditional impression taking.

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33-57506 6mm
33-57508 8mm
33-57510 10mm

bmX®, palatal expander
bmX®, Gaumennahtspreizer


MARPE – Miniscrew Assisted Rapid Palatal Expansion

smartlock®

Video

Quadexpander

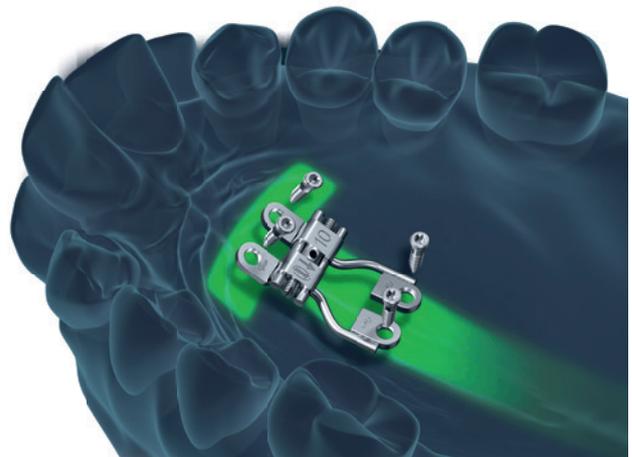BENEFIT® DIRECT Screws | BENEFIT® DIRECT Schrauben

WHY Smartlock BENEFIT DIRECT screws for your perfect MARPE result?

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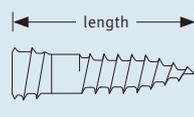
BENEFIT® DIRECT Screw 2.0 mm
BENEFIT® DIRECT Schraube 2,0 mm

2x

ST-33-64209-2	2.0 x 9mm
ST-33-64211-2	2.0 x 11mm

BENEFIT® DIRECT Screws in a package of 2 ea.

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
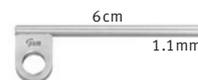
Mouse Trap – for Molar Intrusion | Mausefalle – zur Molarenintrusion

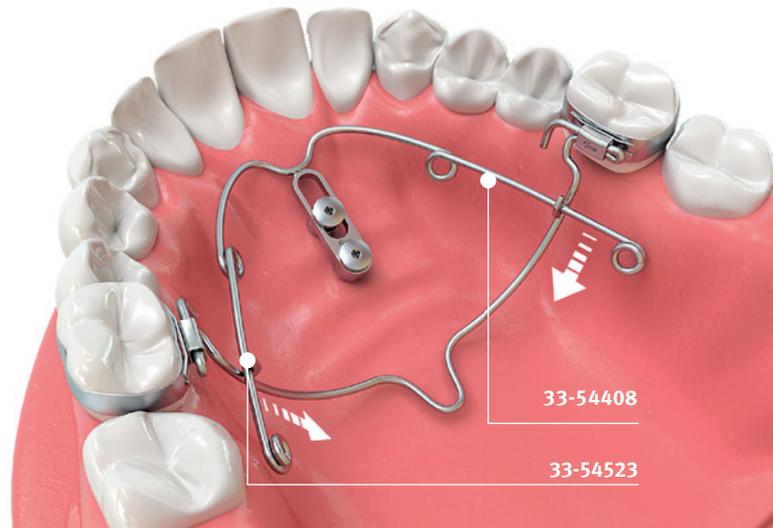consisting of | bestehend aus:

ST-33-54209-2

Page | Seite 8

BENefit® Screw 2.0 x 9mm, 2ea.
BENefit® Screw 2,0 x 9mm, 2Stk.

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Scanbodies | Scanbodies

As a further step into digitalization of the **BENEFIT**® System, these scanbodies were developed in cooperation with Orthodontist Simon Graf (Belp, Switzerland) to ensure visualization of a detailed rendering. The Scanbody Standard is used right after implant insertion. The Scanbody mini can be screwed onto the implants even after several months of treatment.

Developed in cooperation with:
Simon Graf, Orthodontist, Belp, Switzerland

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33-10760

2ea. / 2Stk.

Scanbody PEEK Standard
Scanbody PEEK Standard

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Set of 10 ea.
 10 Stück/Pkg.

Sterile packaging | Sterile Verpackungen

Implants – in single sterile packaging

The outer and inner packaging contains labels indicating UDI codelike article number, LOT number, expiration date, etc. for the patient folder and for full traceability. The sterile screws are contained in a small colour-coded screw carrier which protects the screw and at the same time acts as a pick up tray. This carrier is placed in a clear, sealed blister package. The screw can easily be removed from the carrier by means of the self retaining screwdriver.

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Instruments Instrumente

Drills | Bohrer

The instruments can be used for several **psm** products. In order to find the right instruments for your use, they are indicated in color.

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DENTAL attachment
DENTAL Ansatz

12-12625 1ea. / 1Stk.

Drill, 1.0x33mm WL 15mm, green, for 1.5 / 1.8mm screws
Bohrer, 1,0x33mm AT 15mm, grün, für 1,5 / 1,8mm Schrauben

11-12208 1ea. / 1Stk.

Drill, 1.2x22mm WL 8.5mm, green, for 1.5 / 1.8 mm screws
Bohrer, 1,2x22mm AT 8.5mm, grün, für 1,5 / 1,8 mm Schrauben

10-67513 1ea. / 1Stk.

Drill, 1.4x33mm WL 15mm, red, for 2.0mm screws
Bohrer, 1,4x33mm AT 15mm, rot, für 2,0mm Schrauben

11-18452 1ea. / 1Stk.

Drill, 1.8x33mm WL 15mm, grey, for 2.3mm screws
Bohrer, 1,8x33mm AT 15mm, grau, für 2,3mm Schrauben

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The extension of 16 mm works for all rotary instruments of the **quattro®** and **BENEFIT®**-System with a latch type DENTAL attachment (max. rpm 600).

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
Drill stop | Bohrerstop

The drill stop prevents the clinician from unintentionally drilling too deep. They assure an accurate working length from 7 to 13 mm in 2 mm increments. By sliding the stop sleeve over the appropriate **psm** drill, the working length is adjusted to the desired depth. Removal of the sleeve is made easy by using the buffer.

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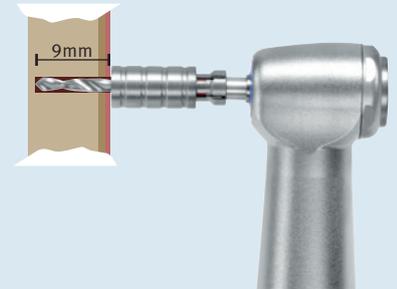33-18307

1 ea. / 1Stk.


Drill-Stop sleeve, 7mm
Bohrer-Stophülse, 7mm

33-18309

1 ea. / 1Stk.

Drill-Stop sleeve, 9mm
Bohrer-Stophülse, 9mm

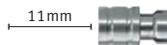33-18311

1 ea. / 1Stk.

Drill-Stop sleeve, 11mm
Bohrer-Stophülse, 11mm

33-18313

1 ea. / 1Stk.

Drill-Stop sleeve, 13mm
Bohrer-Stophülse, 13mm

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33-18288

1ea. / 1Stk.

Laboratory Screwdriver, fixed
Labor Schraubendreher, fest

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33-825-33

1ea. / 1Stk.

BENEFIT® DIRECT TX Screw holder, 28 mm
BENEFIT® DIRECT TX Schraubenhalter, 28 mm

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BENEFIT® Screw | **BENEFIT®** Schraube

BENEFIT® DIRECT Screw | **BENEFIT®** DIRECT Schraube

Battery Screwdriver – iSD900 | Akku Schraubendreher iSD900

Cordless Battery Screwdriver with Torque Calibration System for the controlled fixation of *psm* Miniscrews.

The iSD 900 will guarantee precise torque control due to the unique – easy to use – torque calibration system (10 – 40 Ncm). Time-saving and precise, riskfree insertion and tightening of screws and abutments. The complete set contains a reliable and user friendly handpiece for Dental attachments. Reverse mode is indicated with sound. Standard rechargeable AAA Ni-Mh batteries can be used.

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Fig./Abb.2

Fig./Abb.3

The new Mobilholder is intended to move the **BENEFIT®** mobilizers into the final position prior to tightening the set screw for distalization (fig. 2, 3) or mesialization (fig. 1).

developed in cooperation with ZTM W. Bervoets

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nach ZTM W. Bervoets

Handles | Handgriffe

The manually turned unit is connected to your existing contra-angled handpiece (blue 1:1) that enables its use like an angled screwdriver.

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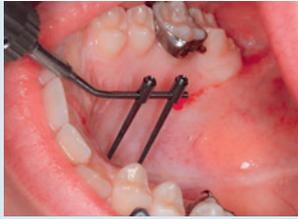
median

paramedian

The soft tissue probes are infinitely adjustable and are fixed by means of the normal **BENEFIT®** instruments (i. e. **33-18288**) with set screws in the desired position.

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The screw location probe acc. to Dr. Schlichter identifies the target points for the median or paramedian placement of the **BENEFIT®** Screws. The precise positioning of the **BENEFIT®** Implants improves the location and parallelism with each other. It also increases stability and simplifies the chairside bending of the plate. Additionally the lab time is reduced due to less bending.

Developed in cooperation with Dr. Kai Schlichter, Dres. Werling & colleagues, Offenbach a. d. Q., Germany.

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
BENEsleeves double, parallel, 2 ea.

BENEsleeves doppel, parallel, 2 Stk.

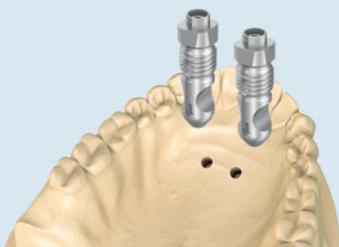The matching drill (**33-10914**) and screwholder with stop (**33-10903**) you will find on page 35.

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
Developed with
Entwickelt mit:
Dr. David J. De Franco
Milano, Italy

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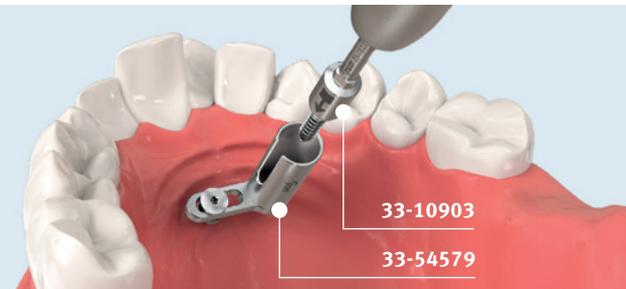The screwholder (**33-10903**) and drill (**33-10914**) can be used with digital motion (page 34), De Franco parallel guide and individually made insertion guides (i.e. TADMAN).

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95-13012 1ea. / 1Stk.

Mouse Trap Intrusion
Mausefalle Intrusion

95-13013 1ea. / 1Stk.

Hybrid Hyrax® Distalizer
Hybrid Hyrax® Distalizer

95-13014 1ea. / 1Stk.

T-Mesial/Distalslider
T-Mesial/Distalslider

95-13001-15 1ea. / 1Stk.

BENEsliider® mesial/distal
BENEsliider® mesial/distal

V-90-994-00 1ea. / 1Stk.

BENEFIT® Model Case for
9 models, empty
BENEFIT® Modellkoffer für
9 Modelle, leer

V-90-994-02 1ea. / 1Stk.

BENEFIT® Model Case for
18 models, empty
BENEFIT® Modellkoffer für
18 Modelle, leer

95-13017-15 1ea. / 1Stk.

BENEsliider® mesial/distal
BENEsliider® mesial/distal

Color | Farbe

Please select the desired **BENEFIT®** model/s and choose the color from the QR code file.

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Catalogues | Kataloge

More innovative **psm** products are presented in our other catalogues. Please request it from your local **psm** distribution partner or download from our website www.psm-medical.com.

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KAT-034

smartlock®

The new **psm – smartlock®** system is largely based on the successful and proven **BENEFIT®** system. **smartlock®** focuses on the essentials: Unique comfort as well as safe and easy handling that doctors and patients can rely on. The trend-setting concept, combined with the best materials and the most modern manufacturing processes, guarantees a treatment result at the highest level.

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
KAT-003

quattro®

The **quattro®** Mini Anchor System is known as one of the first and most successful orthodontic Mini-Implant Systems worldwide. Based on this system the **BENEFIT®** screws were developed. This enables the user to utilize the entire scope of orthodontic implantation with unmatched possibilities and compatible instrumentation.

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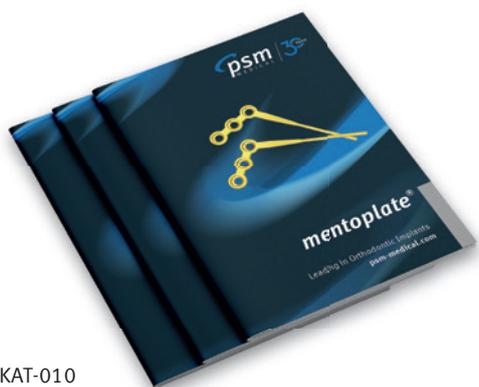KAT-010

mentoplate®

The **mentoplate®** system provides a perfect solution for enhanced stability with special orthodontic bone anchored plates for class III correction. The System includes Titanium plates which are fixed to the bony surface with selfdrilling screws.

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KAT-030

EASY DRIVER®
THE GUIDING SYSTEM FOR TADS

Digital skeletal anchorage | Digitale skelettale Verankerung

A new and innovative method for the secure insertion of the **BENEFIT®** miniscrews. Implantation and appliance placement in one session.

Visit www.uniontech.it/easydriver for more information.

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TADMAN®

The complete digital workflow realized with no efforts for the orthodontist. Based on the **BENEFIT®** System.

Visit www.tadman.de for more information.

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Specifications and design are subject to change without prior notice. All pictures, photos and product descriptions are the property of **psm** Medical GmbH. Utilization and copies by third parties have to be authorized by **psm** Medical GmbH. All rights reserved! All trademarks and brands mentioned in this brochure are subject to unrestricted regulations of the trademark laws valid and right of ownership of the registered owner.

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KAT-002-DE | 09.2023

BENEFIT[®]

Leading in Orthodontic Implants
psm-medical.com

Distribution | Vertrieb:

Manufacturer | Hersteller:

PSM Medical GmbH

Gewerbestr. 10
78594 Gunningen, Germany
Phone: +49 7424 97515-0
E-Mail: info@psm-medical.com
www.psm-medical.com

PSM North America, Inc.

80900 Weiskopf
La Quinta, CA 92253, USA
Toll-free: 800-733-1622
E-Mail: service@psm-na.us
www.psm-na.us