


JV BiTurbo2™

for rapid correction of deep overbite

Ref. A00402

Pack of 10 units

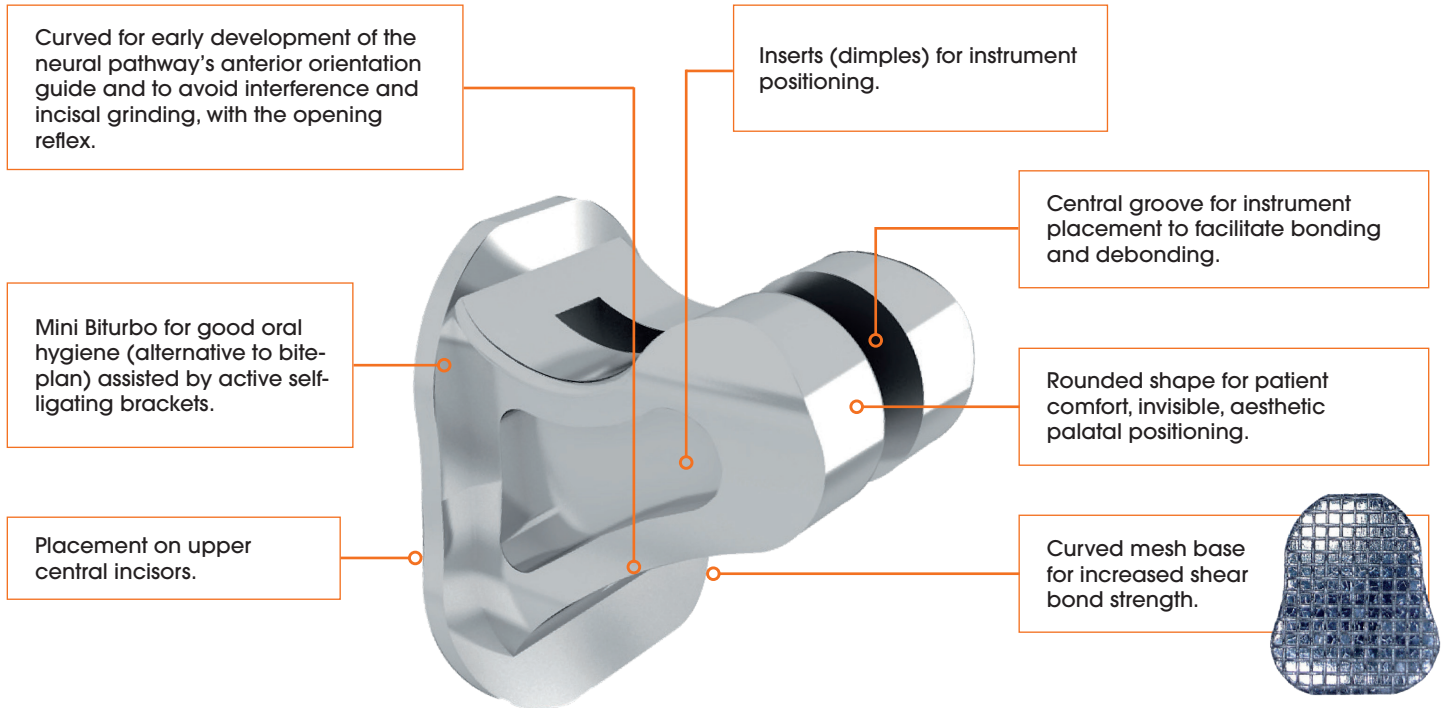

- For rapid correction of severe deep and skeletal supraclussions.
- Reduces hyperactivity and muscle parafunction.
- To be bonded to the palatal surfaces of the upper central incisors to allow posterior disocclusion and avoid interference with the lower incisors.
- Prevents anterior interference.
- Allows physiological eruption of molars and premolars without the effects of molar intrusion.
- Posterior teeth come into contact in 3 to 4 weeks thanks to elastic bands.
- Designed with a central slot to facilitate positioning with a periodontal probe and easy removal by compressing with pliers.
- Made of stainless steel to prevent ingestion of toxic substances.

